



Axial (+3) Image



Axial (+3) Ground Truth

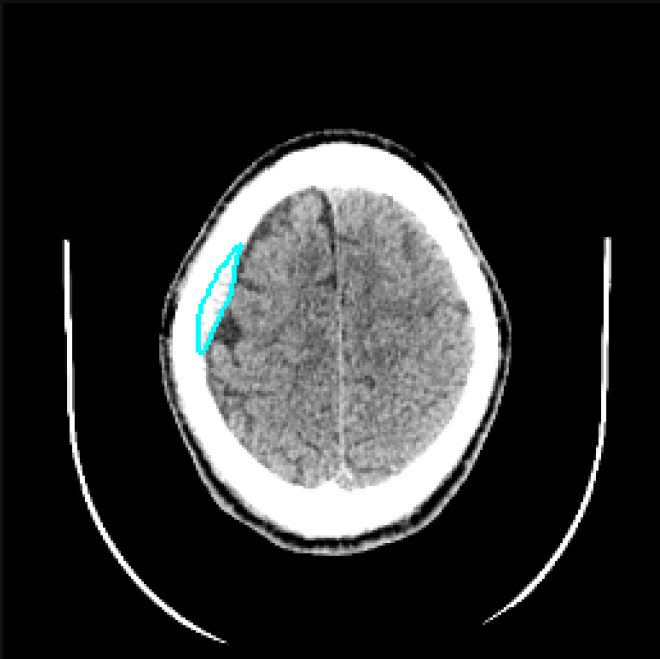

Axial (+3) Eval Composite

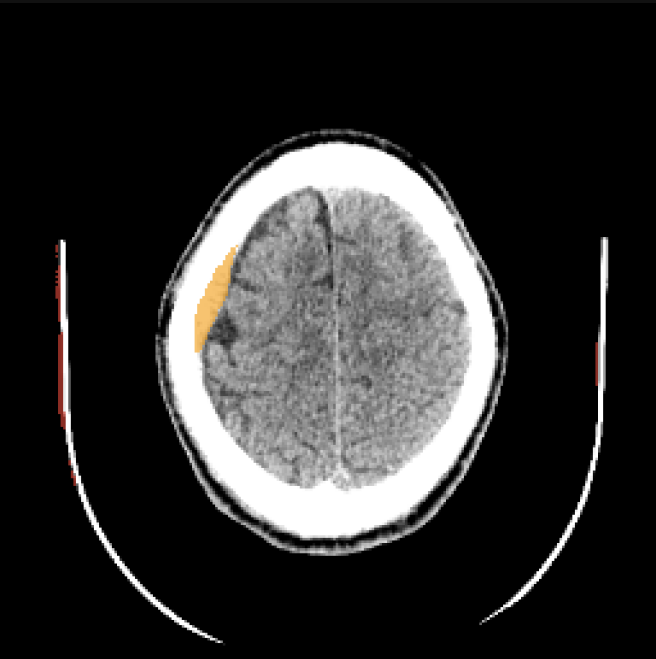

Patient: CQ500_0006
Dice: 0.0000 | IoU: 0.0000
Prec: 0.0000 | Rec: 0.0000 | Spec: 0.9987
TP:0 FP:5480 FN:6075 TN:4182749

GT Vol: 26.25 mL | Pred Vol: 23.68 mL
RVD: -9.8%
HD95: 130.81 mm | ASSD: 62.27 mm


Lesions GT/Pred: 4/9
Lesion Sens/Prec: 0.00/0.00

Optimal Thresh: 0.30 (Dice 0.0000)


Diagnosis:
Severe under-segmentation; missed significant hemorrhage volume.
Multi-focal Case: missed 4 of 4 distinct hemorrhage foci entirely, despite
distant false positive overlap on the detected ones. Boundary distance is
large (HD95=130.8mm), indicating localized regions of substantial
disagreement even where bulk overlap looks fine.

Legend:
True Positive
False Positive
False Negative
Ground Truth


Axial (Center) Image

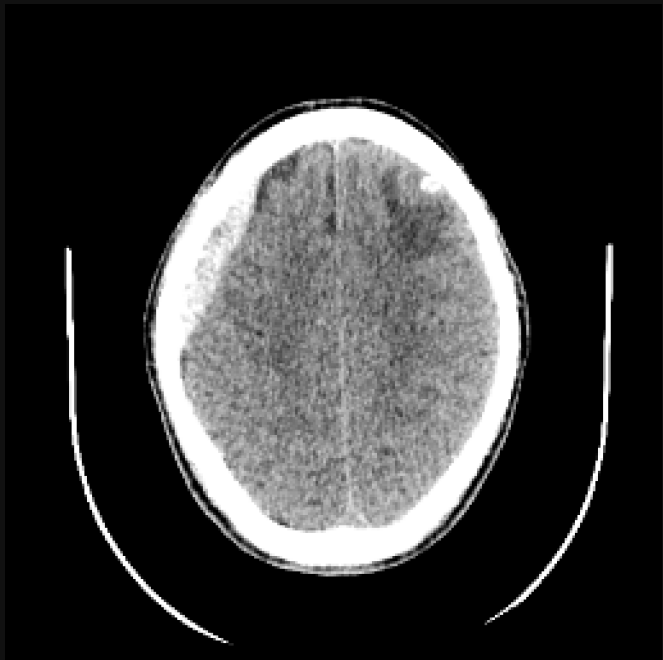


Axial (Center) Ground Truth

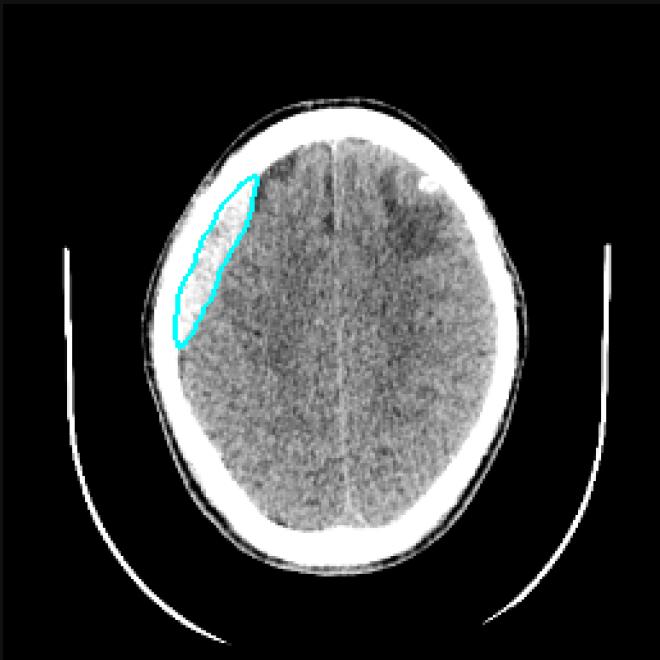


Axial (Center) Eval Composite

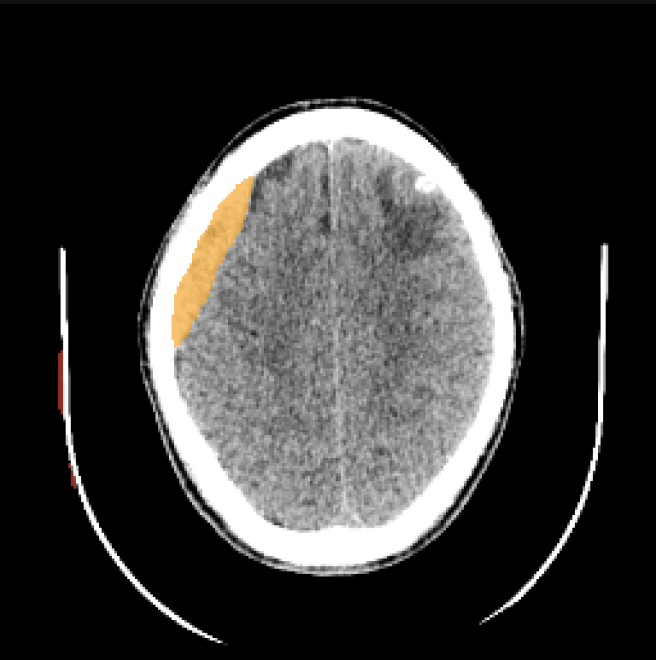

Axial (-3) Image

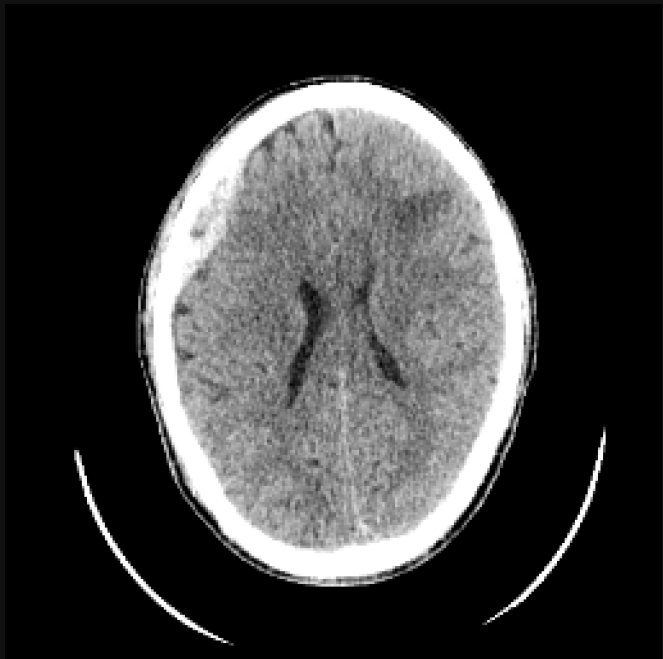

Axial (-3) Ground Truth

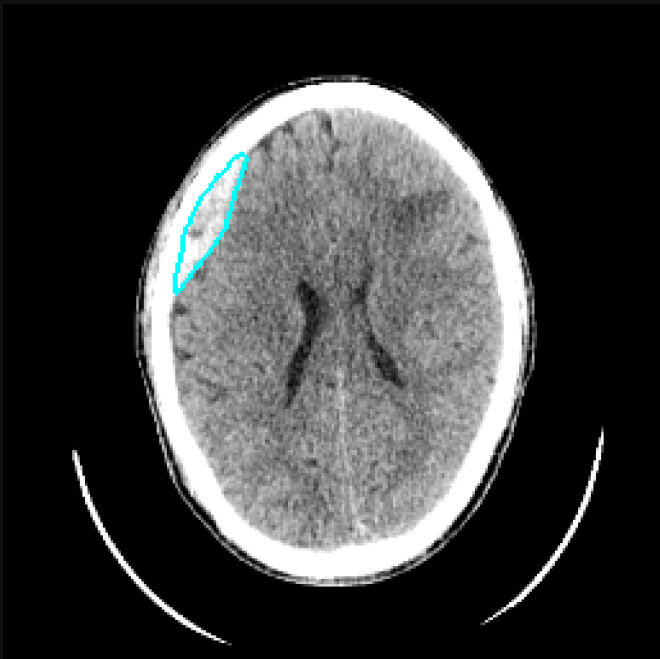

Axial (-3) Eval Composite

