

Axial (+3) Image



Axial (+3) Ground Truth



Axial (+3) Eval Composite



Patient: CQ500_0019
Dice: 0.1312 | IoU: 0.0702
Prec: 0.2803 | Rec: 0.0856 | Spec: 0.9995
TP:794 FP:2039 FN:8479 TN:4182992

GT Vol: 19.37 mL | Pred Vol: 5.92 mL
RVD: -69.4%
HD95: 93.89 mm | ASSD: 61.62 mm

Lesions GT/Pred: 3/4
Lesion Sens/Prec: 0.33/0.25

Optimal Thresh: 0.30 (Dice 0.1466)


Diagnosis:
Severe under-segmentation; missed significant hemorrhage volume.
Multi-focal case: missed 2 of 3 distinct hemorrhage foci entirely, despite
contrastive overlap on the detected ones. Boundary distance is
large (HD95=93.9mm), indicating localized regions of substantial
disagreement even where bulk overlap looks fine.

Legend:
Green: True Positive
Red: False Positive
Blue: False Negative
Red: Ground Truth


Axial (Center) Image



Axial (Center) Ground Truth

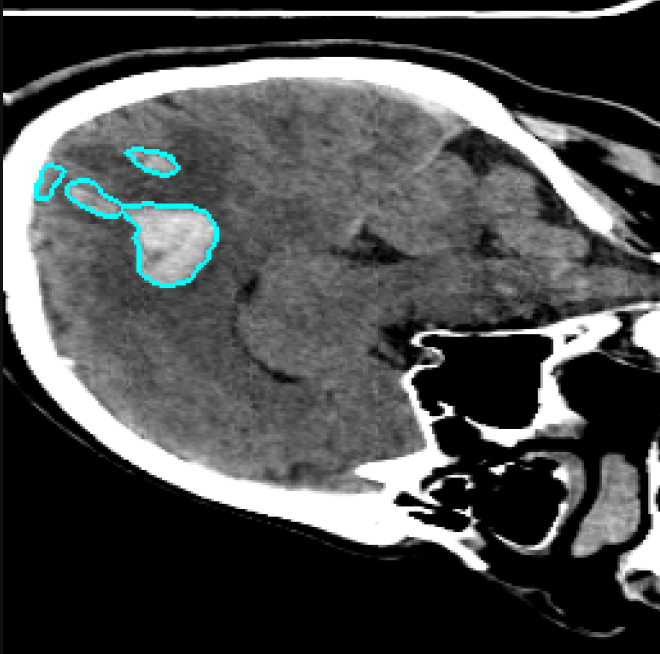

Axial (Center) Eval Composite

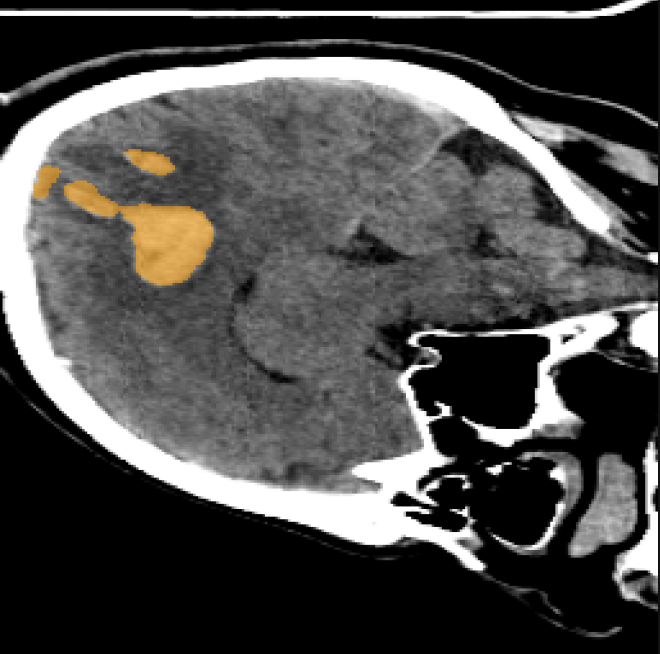

Axial (-3) Image



Axial (-3) Ground Truth



Axial (-3) Eval Composite

