

Axial (+3) Image



Axial (+3) Ground Truth



Axial (+3) Eval Composite



Patient: CQ500_0027
Dice: 0.0000 | IoU: 0.0000
Prec: 0.0000 | Rec: 0.0000 | Spec: 0.9991
TP:0 FP:3609 FN:1385 TN:4189310

GT Vol: 6.60 mL | Pred Vol: 17.21 mL
RVD: +160.6%
HD95: 111.31 mm | ASSD: 87.24 mm

Lesions GT/Pred: 4/5
Lesion Sens/Prec: 0.00/0.00

Optimal Thresh: 0.30 (Dice 0.0000)


Diagnosis:
Severe over-segmentation; spurious false-positive clusters present.
Multi-focal case: missed 4 of 4 distinct hemorrhage foci entirely, despite
distant false-positive overlap on the detected ones. Boundary distance is
large (HD95=111.3mm), indicating localized regions of substantial
disagreement even where bulk overlap looks fine.

Legend:
True Positive
False Positive
False Negative
Ground Truth


Axial (Center) Image



Axial (Center) Ground Truth

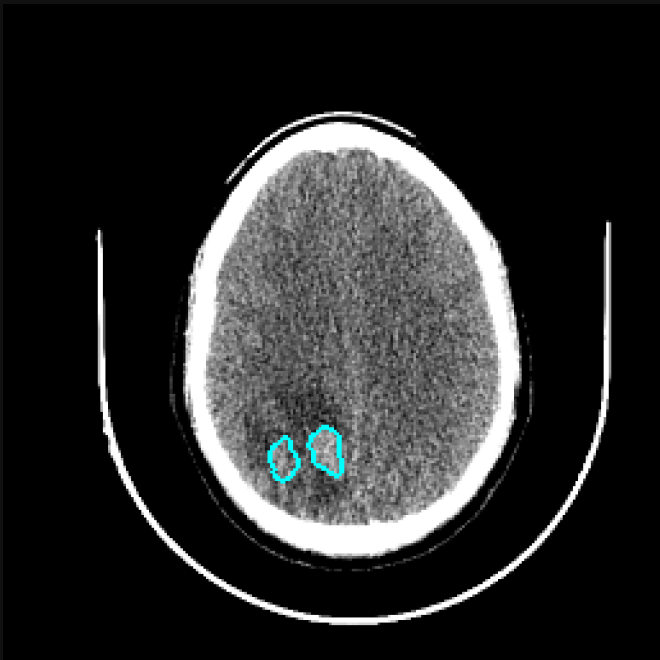


Axial (Center) Eval Composite

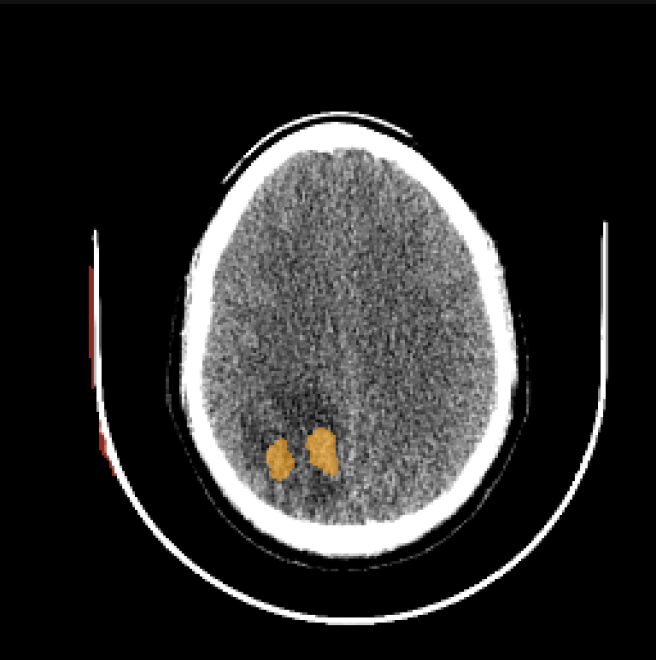

Axial (-3) Image



Axial (-3) Ground Truth

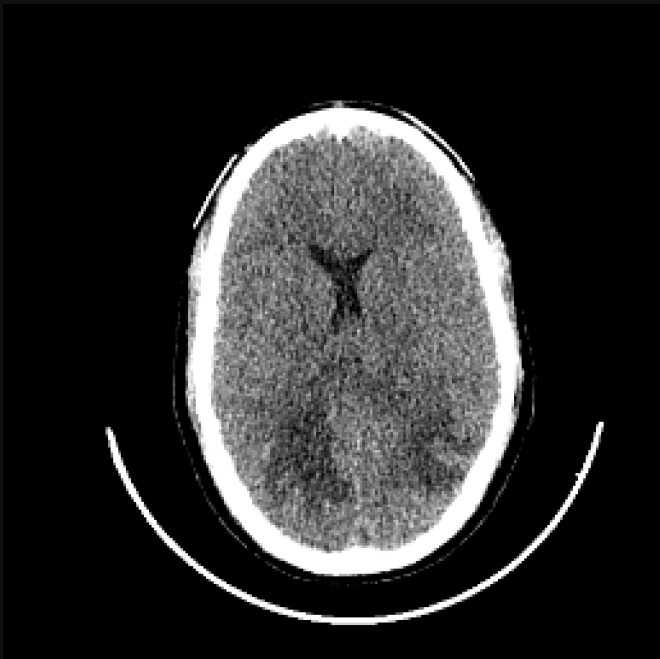

Axial (-3) Eval Composite

